


## PR 39 Pitched-Blade Impeller

P/N: 509-39000-10

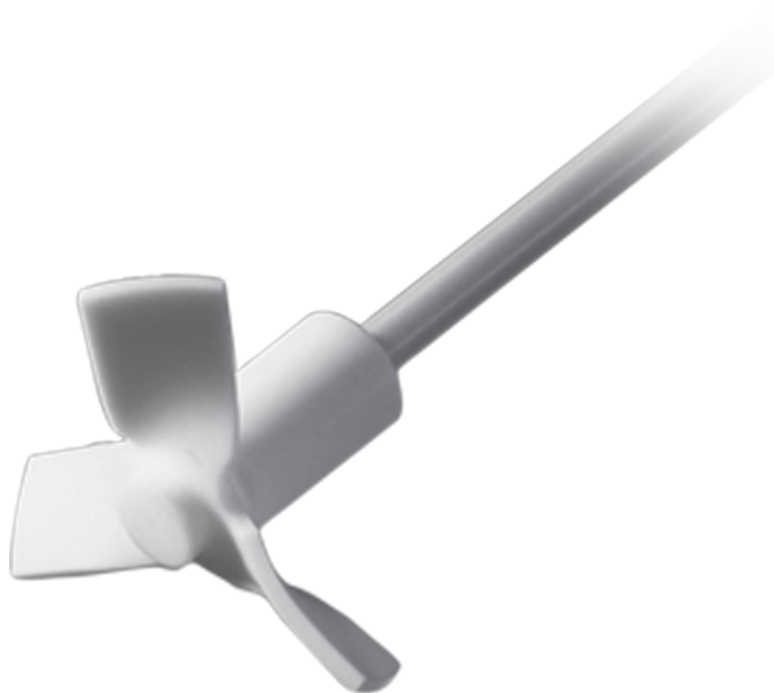

### Propeller-Type Impeller

- These impellers are recommended for applications which require average or high speed
- For mixing tasks with medium or high viscosity
- Excellent mixing properties for homogenization and suspensions
- These models create an axial flow

## PR 39 Pitched-Blade Impeller - Technical Data

|                     |         |
|---------------------|---------|
| Blade size          | 75 mm   |
| Material            | PTFE    |
| Length              | 350 mm  |
| Shaft diameter      | 8 mm    |
| Max. rotation speed | 800 rpm |